Recombinant SARS-COV-2 Spike S1 (N501Y, K417N, E484K) Protein (His Tag)



www.bioss.com.cn 400-901-9800 sales@bioss.com.cn support@bioss.com.cn


Cat. No. bs-46017P

Description	
Protein Sequence	Recombinant SARS-COV-2 Spike S1 (N501Y, K417N, E484K) protein is expressed from mammalian with a His tag at the C-terminal. It contains Val16-Arg685 (N501Y, K417N, E484K).
Source	Mammalian Expression System
Accession	YP_009724390.1
Mol wt	The protein has a predicted MW of 75.9 kDa. Due to glycosylation, the protein migrates to 90-130 kDa based on Bis-Tris PAGE result.
Endotoxin	Less than 1EU per ug by the LAL method.
Purity	> 95% as determined by Bis-Tris PAGE
	>95%as determined by SEC-HPLC
Activity assay	Not tested.
Formulation and Storage	
Formulation	Lyophilized from 0.22um filtered solution in 20mM PB (pH 7.4). Normally 5% trehalose is added as protectant before lyophilization.
Reconstitution	Centrifuge tubes before opening. Reconstituting to a concentration more than 100 µg/ml is recommended. Dissolve the lyophilized protein in distilled water.
Storage	The product should be stored at -70°C or -20°C.
Background	
	SARS-CoV-2 exploits angiotensin-converting enzyme 2 (ACE2) as a receptor to invade cells. It has been reported that the UK and South African strains may have higher transmission capabilities,


SARS-CoV-2 exploits angiotensin-converting enzyme 2 (ACE2) as a receptor to invade cells. It has been reported that the UK and South African strains may have higher transmission capabilities, eventually in part due to amino acid substitutions on the SARS-CoV-2 Spike protein. The results of a study show the N501Y replacement in this region of the interface (present in both the UK and South African strains) should be favorable for the interaction with ACE2, while the K417N and E484K substitutions (South African strain) would seem neutral or even unfavorable.

Assay Data

Tris-Bis PAGE

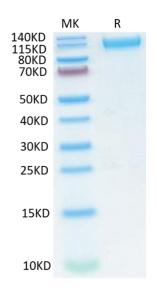

HPLC Data

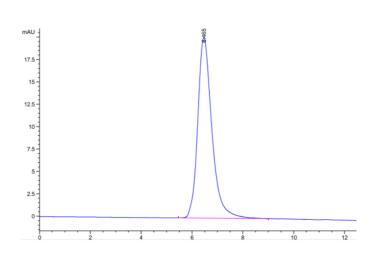

Recombinant SARS-CoV-2 Spike S1 (N501Y, K417N, E484K) on Tris-Bis PAGE under reduced condition. The purity is greater than 95%.

The purity of 2019-nCoV Spike S1 is greater than 95% as determined by SEC-HPLC.